


Emu-Graph 5 Emulator

CI 141

The Emu-Graph 5 chemical emulator strips are designed for use in pre-vacuum steam sterilizers operating at either 121°C for 15 minutes or longer and 134°C for 5.5 minutes or longer. In the presence of quality steam, the Emu-Graph 5 shows the attainment of proper sterilization conditions.

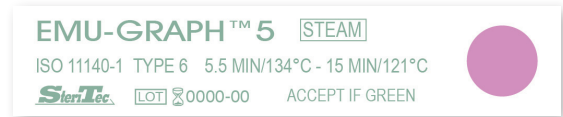

Key Features

- 5.5 minutes at 134°C or 15 minutes at 121°C
- Sterilant specific - must have steam to change
- Lead and toxic heavy metal free

Performance Characteristics

The Emu-Graph 5 indicator changes color when exposed to steam sterilization conditions.

Physical Characteristics

Dimensions: L 89mm x W 18mm

Color Change

Purple to green

Stated Values (as determined in a steam resistometer)

- 6.5 min at 134°C
- 16 min at 121°C

Storage

The storage conditions are printed on the product packaging. Store in the sealed bag provided

Shelf Life

24 months from date of manufacture. The expiry date is printed on the product packaging

Standards Classification

ISO 11140-1, Type 6

SteriTec Part No.

CI 141, Emu-Graph 5 Emulator

Packaging

20 bags of 250