



# Getinge Clean Cleaning Chemistries

The winning formula in the battle against  
mother and human nature

GETINGE 



# The easy choice no matter how hard the water

Mother Nature is not always kind to Central Sterile Department managers. Cleaning and disinfection regimens often are significantly impacted by the scale-producing and corrosive effects of dissolved minerals.

Additionally, even the most stringent processes for washing and disinfection can be a compliance challenge when established procedures for dosing and detergent selection are not followed.

That is why leading institutions in the United States and around the world place their trust in Getinge Clean cleaning chemistries. Getinge's products help ensure their CS staff has every possible advantage to consistently comply with regulatory requirements, while meeting the increasing expectations of hospital administration and a well-informed public.

## **With so much at stake, why take the risk?**

Chances are, over the upcoming years, your institution will invest millions of dollars in advanced CS technology to meet an increasing volume of OR activity and stringent patient safety directives. Getinge Clean products safeguard your investment by helping eliminate the clinical, operational, and financial risks associated with less-than-optimal washing and disinfection.

## **More than 20 years of proven performance**

Getinge's pioneering detergent formulations demonstrate superior cleaning efficacy, even under challenging hard water conditions, while exhibiting excellent corrosion-inhibiting properties. Additionally, Getinge Clean detergents are specially formulated to mitigate the effects of unintended dosage variability.

## **Added confidence comes standard with every container**

Unlike the vast majority of competitive products, Getinge Clean detergents, lubricants, rinse aids, and equipment maintenance cleaners/descalers are produced and packaged in an FDA regulated, medical device facility. This helps ensure consistent quality and performance from batch to batch.

Furthermore, color-coded packaging labels provide immediate visual verification of the product type selected for use. Getinge's color-coded system also facilitates logical product substitutions if a particular chemistry is unavailable. Any other product with the same color-coded label can be substituted with minimal or no changes to equipment cycles.

Packaging also incorporates easy-to-associate symbols that help eliminate uncertainty among experienced and new staff members alike. All containers feature tamper-evident seals in a wide range of sizes for manual dosing and distribution systems.


# Superior formulations for every phase

## Pre-treatment

### Getinge Clean Renuzyme Foam Spray

Foaming, dual-enzyme pretreatment effectively breaks down protein matter from instruments, utensils, and scientific apparatus at point of use. It remains effective even after 96 hours.

Order Number	Description
61301604584	Renuzyme Foam Spray, 12 x 1Qt
61301606712	Renuzyme Foam Spray, 5G

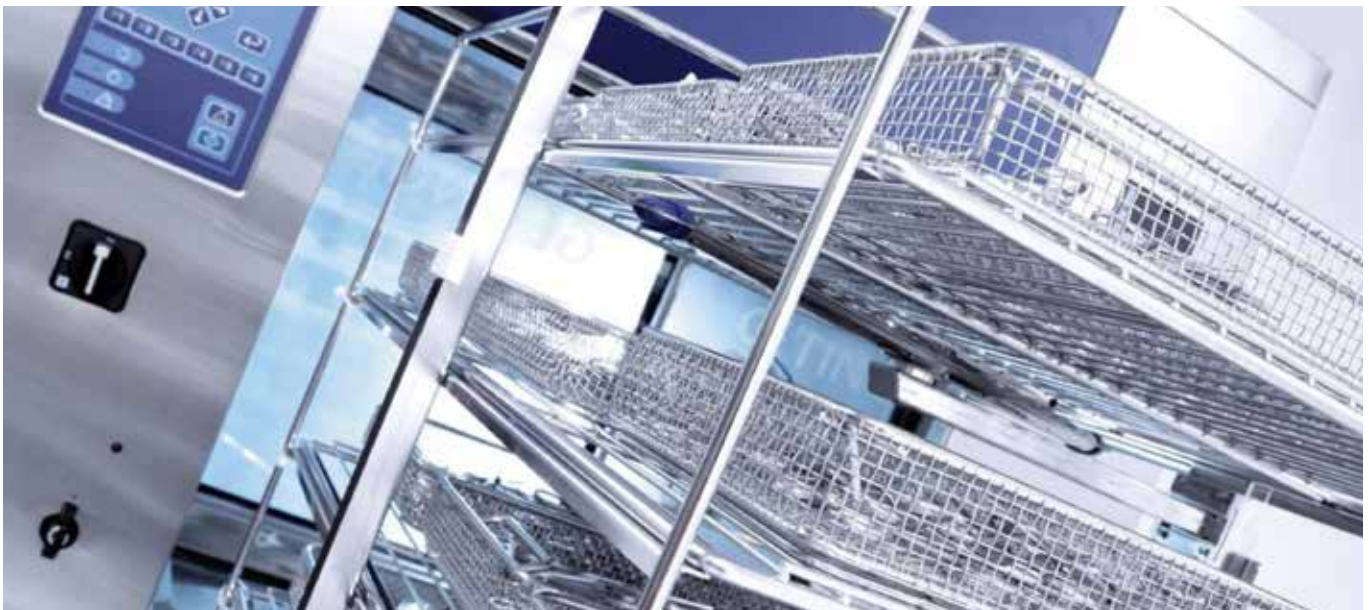

# Manual Cleaning

## Getinge Clean Renuzyme Ultra Enzyme Detergent

Combines advanced protease and lipase enzymes that rapidly digest a broad spectrum of protein and lipid-based soils. Superior wetting action cleans even the most intricate devices. Suitable for sink and machine use.

Order Number	Description
61301606716	Renuzyme Ultra, 4 x 1G
61301606791	Renuzyme Ultra, 2.5G
61301606717	Renuzyme Ultra, 5G
61301606720	Renuzyme Ultra, 15G

## Getinge Clean Renuzyme Plus Detergent

A dual-enzyme formula which quickly penetrates and breaks down protein and organic material and effectively removes blood, tissue, mucous, and other proteinaceous soils in hard or soft water. Suitable for sink and machine use.

Order Number	Description
61301605269	Renuzyme Plus, 4 x 1G
61301600089	Renuzyme Plus, 5G
61301600090	Renuzyme Plus, 15G
61301602690	Renuzyme Plus, 30G
61301601308	Renuzyme Plus, 1 oz. x 50



## Getinge Clean Manual Detergent High Foam

A non-corrosive, concentrated-liquid formulation of anionic and nonionic detergents for manual prewashing of surgical instruments and equipment. The product quickly emulsifies blood and protein for rapid removal. Highly effective in all water conditions and safe for use on stainless steel, aluminum, gold and other soft metals.

Order Number	Description
61301605271	Manual Detergent High Foam, 4 x 1G

*NOTE: Renuzyme WR and PowerCon Triple Enzyme Detergent Concentrate can be used in the manual cleaning phase.*

# Sonic

## Getinge Clean 7900 Sonic Detergent

An effective combination of powerful detergents, wetting agents, and emulsifiers that provide fast and thorough pollution-free cleaning in sonic washers, washer-disinfectors and instrument washers. Non-flammable and phosphate free.

Order Number	Description
61301605272	7900 Sonic Detergent, 4 x 1G
61301600843	7900 Sonic Detergent, 5G
61301601297	7900 Sonic Detergent, 15G

# Detergent Cleaning Chemistries

## Getinge Clean Renuzyme WR Enzyme Detergent

It quickly and effectively removes blood, tissue, mucous, fats, and other soils from surgical instruments and medical devices. For use in sinks, ultrasonic baths, cart washers, instrument washers, and washer-disinfectors.

Order Number	Description
61301605264	Renuzyme WR, 4 x 1G
61301605265	Renuzyme WR, 5G
61301605266	Renuzyme WR, 15G
61301605267	Renuzyme WR, 30G

## Getinge Clean Neutrawash Plus Liquid Detergent

A proprietary blend of detergents, emulsifiers, penetrants, cleaning agents, and water softeners to clean virtually all metals including aluminum. The non-foaming formula works well over a wide range of water conditions.

Order Number	Description
61301605889	Neutrawash Plus, 4 x 1G
61301605890	Neutrawash Plus, 5G
61301605891	Neutrawash Plus, 15G
61301605892	Neutrawash Plus, 30G

## Getinge Clean Tech Wash III Detergent

An economical, one-step aluminum-safe detergent which requires no additives and can be directly dispensed from most automatic systems. Designed for instrument washer-disinfectors and cart washers.

Order Number	Description
61301605275	Tec Wash III, 4 x 1G
61301667776	Tec Wash III, 5G
61301667777	Tec Wash III, 15G
61301600844	Tec Wash III, 30G

## Getinge Clean Acid Detergent

It quickly and thoroughly dissolves minerals and scales from instruments, utensils, and glassware. For use in mechanical cleaning equipment, cart washers, and washer-disinfectors.

Order Number	Description
61301605278	Acid Detergent, 4 x 1G
61301600042	Acid Detergent, 5G
61301600043	Acid Detergent, 15G
61301602229	Acid Detergent, 30G

## Getinge Clean Alkaline Detergent

A highly effective, clear and colorless detergent for removing protein, blood, fat, grease, and oil from surgical instruments, lab ware, and glassware. For use in mechanical cleaning equipment, cart washers, instrument washers, and washer-disinfectors.

Order Number	Description
61301605276	Alkaline Detergent, 4 x 1G
61301600040	Alkaline Detergent, 5G
61301600041	Alkaline Detergent, 15G
61301602228	Alkaline Detergent, 30G



# Concentrated Cleaning Chemistries

## Getinge Clean PowerCon Aluminum Safe Detergent Concentrate

A powerful, low foaming, non-abrasive, detergent that removes, dissolves, and disperses all types of clinical soils, water scale, and oxide stains. Use on a wide range of surgical instrumentation and aluminum trays/pans.

Order Number	Description
61301606043	PowerCon Aluminum Safe Detergent Concentrate, 2 x 1G
61301606044	PowerCon Aluminum Safe Detergent Concentrate, 2.5G
61301606045	PowerCon Aluminum Safe Detergent Concentrate, 5G
61301606694	PowerCon Aluminum Safe Detergent Concentrate, 15G



## Getinge Clean PowerCon Neutral pH Detergent Concentrate

An effective combination of detergents, water conditioners, and wetting agents which delivers superior soil removal over a wide range of operating conditions. Neutral range pH.

Order Number	Description
61301605560	PowerCon Neutral pH Detergent Concentrate, 2 x 1G
61301605561	PowerCon Neutral pH Detergent Concentrate, 2.5G
61301605802	PowerCon Neutral pH Detergent Concentrate, 5G
61301606579	PowerCon Neutral pH Detergent Concentrate, 15G

## Getinge Clean PowerCon Triple Enzyme Detergent Concentrate

This multi-enzymatic compound quickly and thoroughly removes blood, tissue, mucous, fats, and other soils in mechanical cleaning equipment, cart washers, instrument washers, washer-disinfectors, sinks, and ultrasonic baths.

Order Number	Description
61301605563	PowerCon Triple Enzyme Detergent Concentrate, 2 x 1G
61301605564	PowerCon Triple Enzyme Detergent Concentrate, 2.5G
61301605803	PowerCon Triple Enzyme Detergent Concentrate, 5G
61301606581	PowerCon Triple Enzyme Detergent Concentrate, 15G

## Concentrated Cleaning Chemistries *(continued)*

### Getinge Clean PowerCon Acid Detergent Concentrate

A proprietary combination of organic and inorganic acids with surfactants achieves excellent results over a wide variety of operating conditions. The concentrate provides superior dissolving action against minerals and scale. For use in mechanical cleaning equipment, cart washers, instrument washers and washer-disinfectors.

Order Number	Description
61301605569	PowerCon Acid Detergent Concentrate, 2 x 1G
61301605570	PowerCon Acid Detergent Concentrate, 2.5G
61301606233	PowerCon Acid Detergent Concentrate, 5G

### Getinge Clean PowerCon Alkaline Detergent Concentrate

An economical, proprietary formula of powerful alkaline builders, water conditioners, anti-redepositing agents, and surfactants provides high-performance cleaning and effective removal of protein, blood, fat, grease, and oil from surgical instruments, lab ware and glassware. For use in mechanical cleaning equipment, cart washers, instrument washers, and washer-disinfectors.

Order Number	Description
61301605566	PowerCon Alkaline Detergent Concentrate, 2x1G
61301605567	PowerCon Alkaline Detergent Concentrate, 2.5G
61301605800	PowerCon Alkaline Detergent Concentrate, 5G





## Lubricants and Rinse Aids

### Getinge Clean Instrument Lubricant Plus

Extends instrument life and helps ensure proper functioning of hinges, pins, joints, and box locks. No rinsing or wiping required. Steam and EtO permeable.

Order Number	Description
61301605247	Instrument Lubricant Plus, 4 x 1G
61301606790	Instrument Lubricant Plus, 2.5G
61301605248	Instrument Lubricant Plus, 5G
61301605249	Instrument Lubricant Plus, 15G
61301605250	Instrument Lubricant Plus, 30G

### Getinge Clean Tec Rinse

Injects into the final, fresh-water rinse of automatic cart and utensil washers to reduce drying time and prevent spotting in poor water conditions. Compatible with many chemical sterilizers.

Order Number	Description
61301605280	Tec Rinse, 4 x 1G
61301606792	Tec Rinse, 2.5G
61301664156	Tec Rinse, 5G
61301664157	Tec Rinse, 15G
61301601299	Tec Rinse, 30G



## Disinfectants

### Multi-Quat tb

A highly effective, one-step disinfectant cleaner and deodorant for hard, non-porous environmental surfaces within patient rooms, ORs, ICUs, public restrooms, and other areas where housekeeping is challenged to control the hazard of cross contamination. Combination tuberculocidal, disinfectant, mildewstat, fungicide, virucide, cleaner and deodorizer.

Order Number	Description
61301606713	Multi-Quat tb, 12 x 1Qt

### Tec-Quat 128

An economical, high-performance disinfectant, cleaner, sanitizer, fungicide, mildewstat, virucide, and deodorizer. It significantly reduces the potential for cross-contamination on hard, nonporous surfaces. Non-acidic, non-scratch concentrate cleans, disinfects, and deodorizes in one labor-saving step.

Order Number	Description
61301605283	Tec-Quat 128, 4 x 1G
61301605064	Tec-Quat 128, 5G



# Equipment Maintenance

## Getinge Clean Instrument Brightener

Removes rust, corrosion, and hard-water deposits from stainless steel instruments; loosens box locks and joints; and mitigates minor pitting to extend instrument life. Restores the bright original finish to the instrument surface.

Order Number	Description
61301605067	Instrument Brightener, 12 x 1Qt

## Getinge Clean Stainless Steel Chamber Cleaner

This acid detergent cleans and brightens the chambers and exteriors of stainless steel equipment. It dissolves scale and discoloration from calcium carbonates, magnesium carbonates, and sulfates.

Order Number	Description
61301605284	Stainless Steel Chamber Cleaner, 4 x 1G

## Getinge Clean Stainless Steel Polish

A soy-based formulation that removes dirt, grease, fingerprints, and other soils from stainless steel. For use on mechanical cleaning equipment, cart and equipment washers, washer-disinfectors, and sterilizers.

Order Number	Description
61301606201	Stainless Steel Cleaner/Polisher, 1 x 16 oz.

## Getinge Clean Chamber Cleaner Schedule 7

A multipurpose, concentrated acidic cleaner/descaler combines organic and inorganic acids with a controlled foaming surfactant, inhibitors, rinse aids, and penetrants to remove water scale and mineral deposits.

Order Number	Description
61301605287	Chamber Cleaner Schedule 7, 4 x 1G

## Getinge Clean Chamber Cleaner Schedule 90

A concentrated, cleaner/descaler which removes scale and mineral deposits by dissolving binding carbonates. For use inside of sterilizers, mechanical cleaning equipment, cart washers, and cage washers.

Order Number	Description
61301605286	Chamber Cleaner Schedule 90, 4 x 1G

## Getinge Clean Boiler Descaler

A non-corrosive and non-flammable formula which is fortified with wetting and penetrating agents for the rapid removal of water scale, lime, and rust deposits found on water-cooled, water-heated, or water-operated boilers.

Order Number	Description
61301605285	Descaler, 4 x 1G





**Ordering Information****Phone:** 800 950 9912**Fax:** 800 950 2570**Email:** cs.purchaseorders@getinge.com

Getinge is a registered trademark of Getinge AB or its subsidiaries or affiliates in the United States or other countries.  
Copyright 2019 Getinge AB or its subsidiaries or affiliates. • All rights reserved.

**Sales Office, US** • 45 Barbour Pond Drive • Wayne, NJ 07470

**Sales Office, Canada** • 90 Matheson Blvd West, Suite 300 • Mississauga, Ontario L5R 3R3 • Canada

**[www.getinge.com](http://www.getinge.com)**

MCV00058545 REVB  
EXP11/21