


# Wash-Checks Holder

## WC 102

The Holder for the Wash-Checks and Wash-Checks U is an auxiliary device for routine monitoring of automated wash processes in washer disinfectors and ultrasonic washers. The Wash-Checks Holder should only be used with the Wash-Checks (in the washer-disinfector) or the Wash-Checks U (in the ultrasonic washer). The Wash-Checks or Wash-Checks U is inserted into the Wash-Checks Holder, simulating a contaminated surgical instrument.

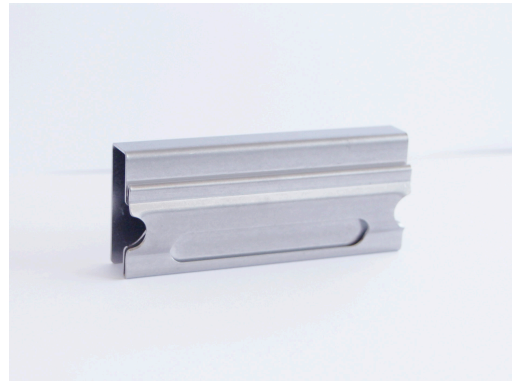

### Key Features

- Reusable holder
- Passivated stainless steel

### Performance Characteristics

The Wash-Checks is inserted into the holder so that half of the test soil is obscured and the other half exposed. The Wash-Checks U is inserted into the holder so that the test soil is fully exposed.

### Physical Characteristics

- **Dimensions:** 89 mm x 38 mm x 15 mm
- **Weight:** 95 g

### Shelf Life

Not Applicable

### SteriTec Part No.

WC 102, Wash-Checks Holder

### Packaging

6 each